

---

# Ultrafiltered Yeast Hydrolysate

---

Catalog No. 96863

1 kg Powder

---

## INTENDED USE

Ultrafiltered Yeast Hydrolysate may be used as a supplement to improve growth and production of a variety of cell cultures.

---

## QUALITY ASSURANCE

All Quality Control test results are reported on a lot specific Certificate of Analysis which is available upon request.

---

## STORAGE INSTRUCTIONS AND STABILITY

Handle using aseptic techniques to avoid contamination.

Store dry at 15-30°C.

Do not use after the assigned expiration date.

---

## SUGGESTED INSTRUCTIONS FOR USE

1. Add 2-6 g/L of hydrolysate supplement to 1 L of culture medium (optimal amount may vary depending on cell line and chosen hydrolysate)
  2. Dissolve completely.
  3. Sterile filter supplemented medium by 0.2  $\mu$  membrane filtration.
-

---

**ALSO AVAILABLE FROM IRVINE SCIENTIFIC**

Optimized Hydrolysate Mixture for CHO Cells Catalog No. 96850

Ultrafiltered Soy Hydrolysate Catalog No. 96857


IS CHO-CD™ Catalog No. 91119

IS CHO-CD4™ Catalog No. 91100

IS CHO-CD and IS CHO-CD4 are chemically-defined media developed for high-density culture of Chinese Hamster Ovary (CHO) cells.

IS MAB-CD™ Catalog No. 91104

IS MAB-CD is a chemically-defined medium formulated for the growth of hybridoma and myeloma cell lines for recombinant monoclonal antibody production.

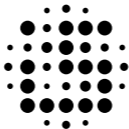

**IrvineScientific®**

2511 Daimler Street, Santa Ana, California 92705-5588

Telephone: 1 949 261 7800 • 1 800 437 5706

Fax: 1 949 261 6522 • [www.irvinesci.com](http://www.irvinesci.com)

PN 40804 Rev. 1

---