

PRODUCT SHEET

IDENTIFICATION OF PRODUCT

NAME OF PRODUCT	PETG EMBRYO STRAW, RIGID
REF. NUMBER	006433, 006578
PACKAGING	BLISTER OF 4 X 5 UNITS
APPLICATION	PACKAGING AND CRYOGENIC STORAGE OF EMBRYOS
COUNTRY OF ORIGIN	FRANCE
MANUFACTURER	CRYO BIO SYSTEM, GROUP IMV TECHNOLOGIES
MEDICAL DEVICE	CE CLASS 2A, STERILE
STERILIZED	YES (STERILIZATION BY IRRADIATION)
SHELF LIFE	3 YEARS

COMPOSITION/INFORMATION ON THE COMPONENTS

Tube = PETG (PolyEthylene Terephthalate), biocompatible, non toxic

Plug = Polyester viscose + polymerizing powder

SPECIFICATIONS

Length = +/- 91 mm for Ref. 006433; +/- 133 mm for Ref. 006578

Diameter (outer)= +/- 1.9 mm

BIOCOMPATIBILITY

Validated on human embryos

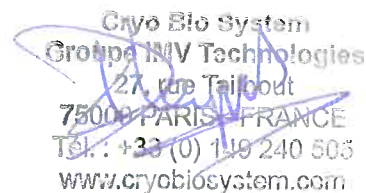
MANIPULATION AND STORAGE/STABILITY

Precaution for Use : The powder of the plug should be wet between 1/3 and 1/2

The storage of the straw in LN2 should be cotton plug at the bottom, ID rod at the top

Do not use if the product or the packaging is damaged as sterility and integrity of the product can no longer be guaranteed

Storage precautions: Preserve away from light, heat and damp


Cryo Bio System
Groupe IMV Technologies
27, rue Tailbout
75009 PARIS - FRANCE
Tél. : +33 (0) 149 240 503
www.cryobiosystem.com

CRYO BIO SYSTEM IS AN ISO 13485:2003 CERTIFIED COMPANY